

Crestpoint ROUNDS™

JANUARY
2019
Volume 3
Issue 1

PUBLISHED FOR ALL OPHTHALMIC AUDIENCES: SURGEONS, YOUNG EYE SURGEONS, NURSES/TECHNICIANS AND PURCHASING TEAMS

IN THIS ISSUE

GEUDER.....	1, 3
Duckworth & Kent.....	2
MANI.....	1
Crestpoint AAO 2018 Recap.....	4

Lamellar and Penetrating Keratoplasty (Corneal Transplantation)

Featuring **GEUDER NANOEDGE Single-Use Trephines and Corneal Punches**

Geuder® NanoEdge
Precision made in Germany

NANOEDGE Single-Use Corneal Punches **NANOEDGE Single-Use Vacuum Corneal Punches**

NANOEDGE Single-Use Cutting Block **NANOEDGE Single-Use Vacuum Trephines** **NANOEDGE Single-Use Trephines**

MANI SAFETY KNIFE

MANI®

Ophthalmic Knife With Safety Cover

SINGLE-USE, DISPOSABLE

Stainless steel
Bevel shape
Non-glare
Crosslapped edge
Ultimate sharpness

AVAILABLE TYPES:

Safety Straight Knife (MST-SK)

Safety Slit Knife (MSL-SK)


Safety MVR Knife (MVR-SK)

Safety Crescent Knife (MCU-SK)

Product shown is
MSL22SK
Slit Angled 2.2mm



MANI Ophthalmic Knife

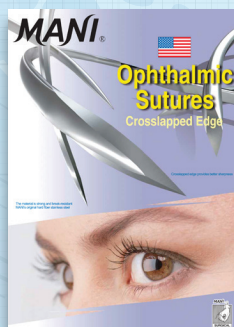


MANI Ophthalmic Sutures




Sterimedix
Main Catalog



Duckworth & Kent
Main Catalog

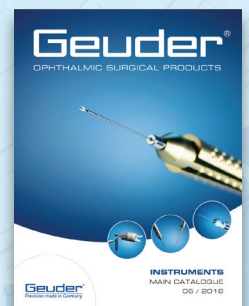

Geuder
Main Catalog

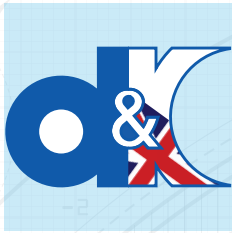

Duckworth & Kent Titanium. Reusable.

Original manufacturer of titanium surgical instruments for ophthalmology and known throughout the world for quality and innovation. D&K titanium instruments are lightweight, non-magnetic, and anodized to reduce glare under the microscope.

Suggested Corneal Graft Set

This set of instruments is compiled to facilitate the surgeon in performing penetrating keratoplasty (corneal transplantation)



1-401 - Castroviejo Corneal Scissors (curved to left)



2-500E - McPherson Tying Forceps



4-115 - DK Retractable Diamond Knife, 1mm Trifacet




9-732 - Thornton Lowest Profile Parallax Free Blade Radial Marker



Visit

www.duckworth-and-kent.com
and navigate to our literature library to view our Suggested Corneal Graft Sets

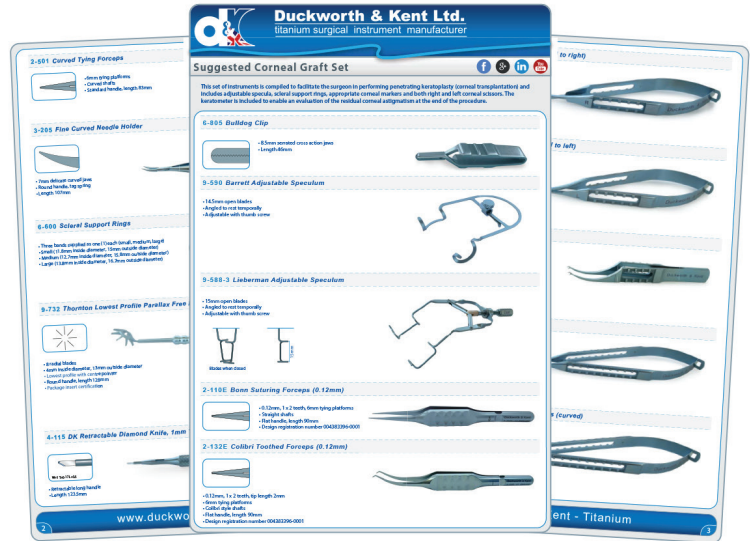

at the leading edge



Years

Duckworth & Kent are proud to celebrate 60 years of being at the leading edge of manufacturing.

1959 - 2019





GEUDER[®] Chromium Plated Stainless Steel. Reusable.

Leading manufacturer of high quality surgical instruments for ophthalmology, made with chromium plated stainless steel. It has always been the maxim of GEUDER to guarantee the superior level of craftsmanship, precision and reliability.

NANOEDGE SINGLE-USE TREPHINES AND CORNEAL PUNCHES NanoEdge



NANOEDGE Single-Use Vacuum Trephines

- 16 radial suture markers, non-slip finger grips
- Seamless blades for precise trephination
- Cross hairs for accurate alignment with the visual axis
- Cutting depth control: 0.0625mm blade progression per quarter turn and 0.25mm per full 360° rotation for precise cuts in various corneal layers
- Wide range of different sizes in small diameter increments for nearly all corneal applications

NANOEDGE Single-Use Trephines

- Seamless blade for precise cuts and a controlled post-operative result
- 0.25mm diameter increments on sizes between 6.0mm to 9.0mm
- Wide range of different sizes in small diameter increments for nearly all corneal applications

NANOEDGE Single-Use Vacuum Corneal Punches

- Vacuum Corneal Punches
- Vacuum fixates the graft and allows for a perfect circular cut of the corneoscleral donor button
- Seating Ring
- Pressing the cornea down to the cutting block (before cutting) supports a non-slip preparation of the graft



NANOEDGE Single-Use Corneal Punches

- Seamless blade for precise punches
- Maximum contrast between graft and white cutting block
- Guide posts support precise tissue cuts without slipping or tilting
- Non-slip positioning of the graft through seating ring and application of vacuum (with G-560xxx)

NANOEDGE Single-Use Cutting Block


- Teflon
- Centric hole facilitates positioning and fixation of the graft during preparation
- Maximum contrast between graft and white cutting block allows precise positioning of the corneoscleral graft to achieve best trephination results



[Click For GEUDER NANOEDGE Brochure](#)

YAMANE DOUBLE-NEEDLE STABILIZER

Globe Fixation Ring and Needle Guide for Flanged Intrascleral IOL Fixation


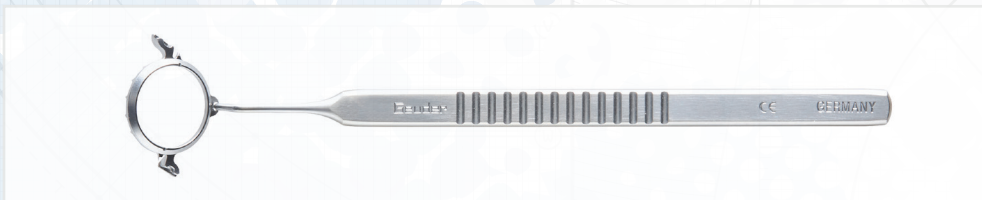


G-31496 YAMANE Double-Needle Stabilizer 30G


G-31497 YAMANE Double-Needle Stabilizer 27G

- Fixation ring for intrascleral IOL fixation with double needle technique
- Overall length 110mm
- The instrument incorporates two landmarks for visualization of the sclerotomy sites, located at 2mm from limbus and 180° apart. Eliminating the need of axis markers, calipers, or ink.

[Click For Brochure](#)

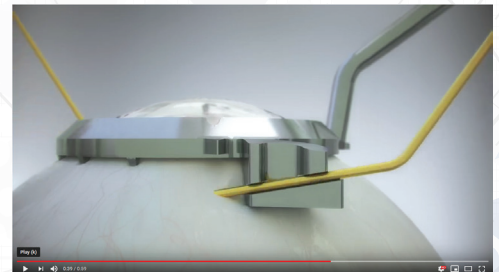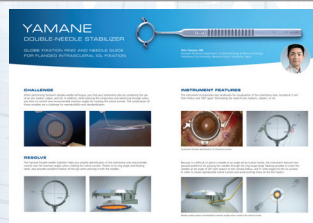

[Click For Video](#)

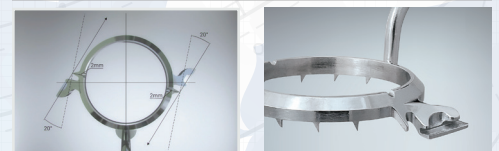

AAO 2018 Recap

Thank you to all who visited Crestpoint Ophthalmics at the AAO 2018 Chicago, October 27-30 – McCormick Place. Our team arrived on Thursday and Friday for setup, and was ready when the first guests arrived on Saturday morning to start a very busy show. We had a great turnout of customers, colleagues and friends from the USA, and abroad.



Contact Us

Crestpoint Ophthalmics is dedicated to being the foremost, industry leading, single-source provider of the world's highest quality ophthalmic products.

Our top priority is superb customer service to you.



Crestpoint Management
Crestpoint Ophthalmics
9100 Watson Road
St. Louis, MO 63126 USA

- Phone: (800) 756-7773
- Phone: (314) 849-7773
- Fax: (314) 849-7776
- E-mail: custsvc@crestpointmgt.com
- Web: www.crestpointmgt.com



GSA/FSS Contract
V797P-4438B
Crestpoint Instruments
SIN A8a & A8b
Effective
5/1/2011 to 4/30/2021

Published by Crestpoint Management. Manufacturer trademarks, part data, and copyrighted materials are property of their respective owners and are used with permission. D&K, Geuder, MANI product ranges in this publication are exclusive in the USA / NA.

Copyright 2019 Crestpoint Management / Crestpoint Ophthalmics
All rights reserved. All information in this publication is subject to change without notice.

ZCC483NL (01/2019)